

MedSafe Extension Sets

MedSafe extension sets are designed with the smallest patients in mind. With neonatal patients, drug incompatibility becomes a considerable liability. Reduce the risk by using extension sets with the least amount of dead space to ensure that your patient is truly receiving the medication. Sets with larger dead space will allow the medication to become trapped in the extension set rather than being infused to the patient.

PRODUCT FEATURES & BENEFITS:

MINIMAL DEAD SPACE: Our quadfurcated connector is designed to minimize potential drug incompatibility in the fluid path with the common fluid juncture located at the male luer (multiple small volume medications delivered in a common leg increases the risk for incompatibility problems and pooling of the medication in the line).


REDUCED CHANCE OF INFECTIONS: Our set is complete and reduces the risk of infection by minimizing the need of additional sets thereby reducing touch contamination.

LOW PRIMING VOLUMES: Our sets have very low priming volumes and short lengths.

LIGHTWEIGHT: Our connector is nearly weightless with an easy rotating luer lock. This will reduce the chance for kinking and pulling.

COLOR CODED: Slide clamps are different colors for easy line identification.

Non-DEHP, Latex free.

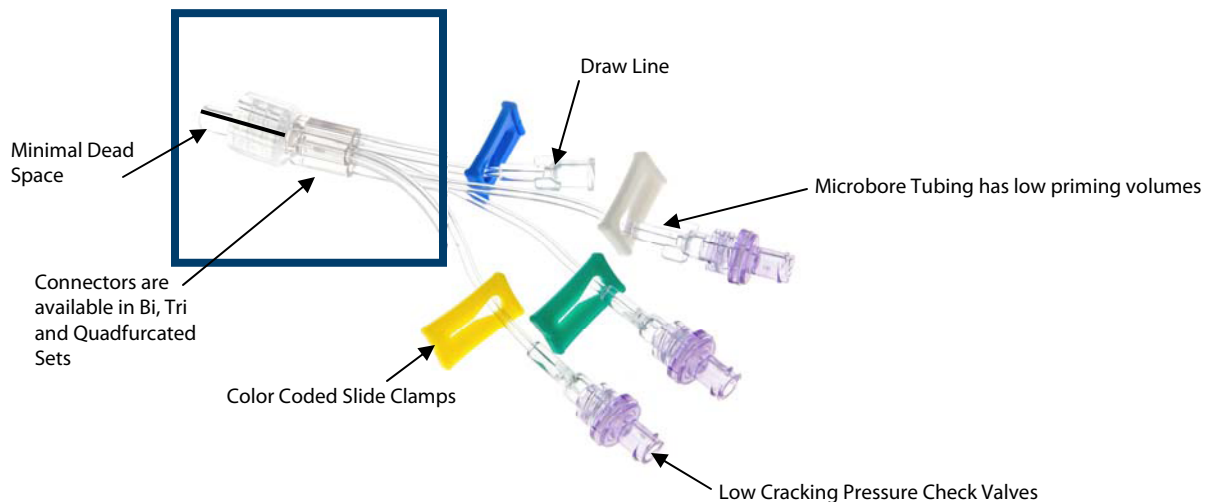

Acacia®

Ph 714.257.0470

Fx 714.257.0513

www.acaciainc.com

IMA
Independent Medical Associates

Products distributed by Independent Medical Associates: www.I-MA.com · info@I-MA.com · p: 888.548.4462 · f: 888-548-1462