


# PILOT Round Tip™ Needle

Precise, safe needle placement is critical to successful transforaminal epidural steroid injections. PILOT's non-cutting, round-tip needle design provides numerous benefits that can help alleviate problems associated with both sharp and blunt tip needles.

- PILOT's round-tip needle is designed to displace tissue along the pathway to the targeted nerve— minimizing the chance intraneural or intravascular injury.
- The round-tip needle allows for smooth placement closer to the targeted nerve or nerve root.
- Medication port is located at the tip of the needle for more effective, precise drug delivery.
- Round tip enhances tactile feel.
- Stylet can be quickly and easily removed with a small twist. Cap is larger for better handling.
- Available in color-coded 22 and 25 gauge sizes and includes introducer needle.



*Illustrated detailed view of Pilot's round tip.*

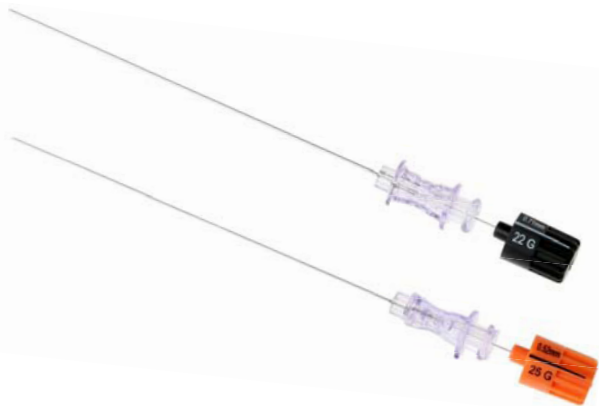

Catalog Number	Gauge	Length	Hub Color	Qty/Box
NB-2215*	22	1.5"	Black	25
NB-2240	22	4.0"	Black	25
NB-2260	22	6.0"	Black	25
NB-2515*	25	1.5"	Orange	25
NB-2540	25	4.0"	Orange	25
NB-2550	25	5.0"	Orange	25

*Sterile Single-Use only. Refer to Directions for Use for complete instructions  
\* Catalog numbers NB-2215 and NB-2515 are Build to Order items.  
Please allow for 8-10 weeks.*