

Critical Care Products

Bronchitrac® "L" Set

TRAKE-fit™

COMFIT™

NEO-fit™

**Esophageal Balloon
Catheter Set**

**INCA® Nasal
Cannulae Set**

NEO-prep



ACKRAD LABORATORIES, INC.
a CooperSurgical Company

Bronchitrac page 1 TRAKE-fit page 2 COMFIT page 3 NEO-fit page 4
Esophageal Balloon Set page 5 INCA page 6 NEO- prep page 6

Bronchitrac® "L" Set



Automatic entry into the left main-stem bronchus

The only suction catheter specifically designed for left-bronchus access in an easy-to-use tray.

Simply advance the Bronchitrac "L" through a properly placed endotracheal or tracheostomy tube. There is no need to twist the catheter or patient. The dual plane, patented curvature of the catheter assures consistent entry into the left main-stem bronchus. The time savings and reduced patient trauma provided by the Bronchitrac "L" makes it the obvious choice for difficult left-bronchial suctionings.

- Ideal for suctioning of post-cardiac surgical, neuro and critical care patients
- Can be used for deep-lobe suctioning
- Completely latex-free product

These catheters are not recommended for right bronchial suctioning.

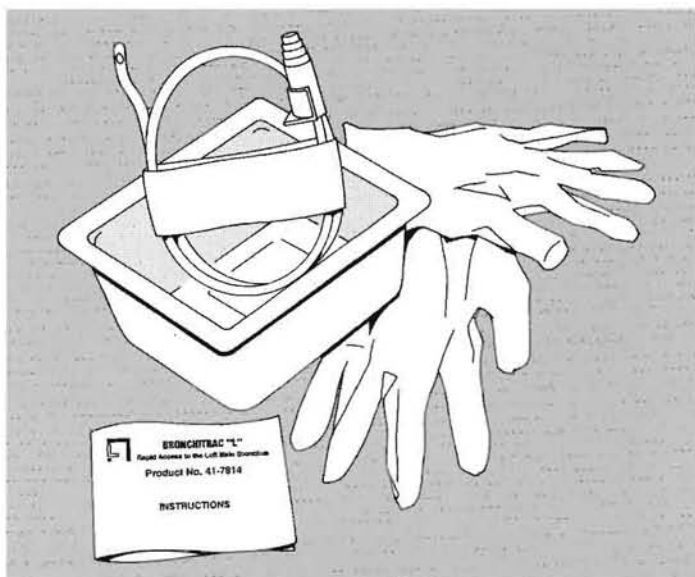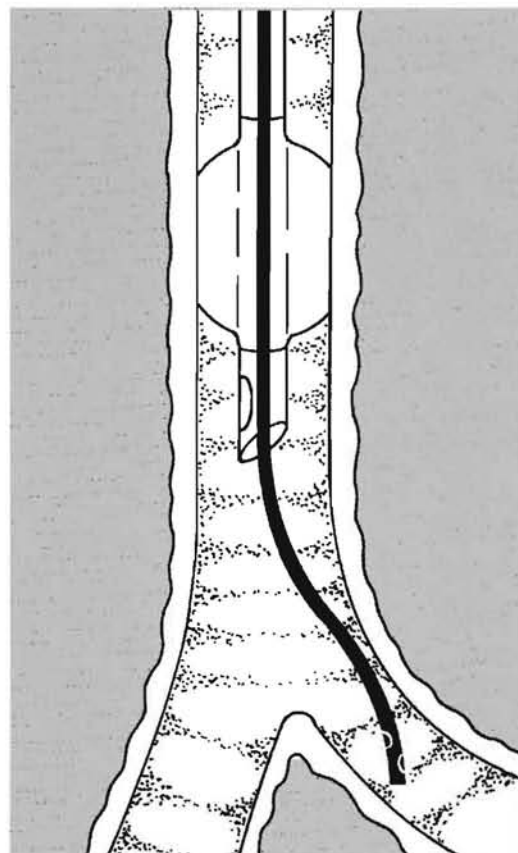

Ordering Information

Bronchitrac "L" Set

Product No. 41-7814

Includes 21 inch long, 14F left-bronchial suction catheter with end hole and two side holes at tip, slide-cover proximal connector, and two vinyl gloves. Contained in clear tray.



Packaging

Sterile, 20 per box, 80 per case (4 boxes).

TRAKE-fit™

Tracheostomy tube holder for maximum patient comfort

- Available in single or double collar sets
- Installs in seconds, remains secure
- Maximum comfort with cotton-lined neckband
- No adhesive, no narrow tape
- Completely adjustable
- Narrow tabs available for pediatric and small adult use
- Third-party reimbursable

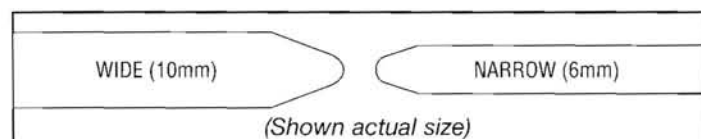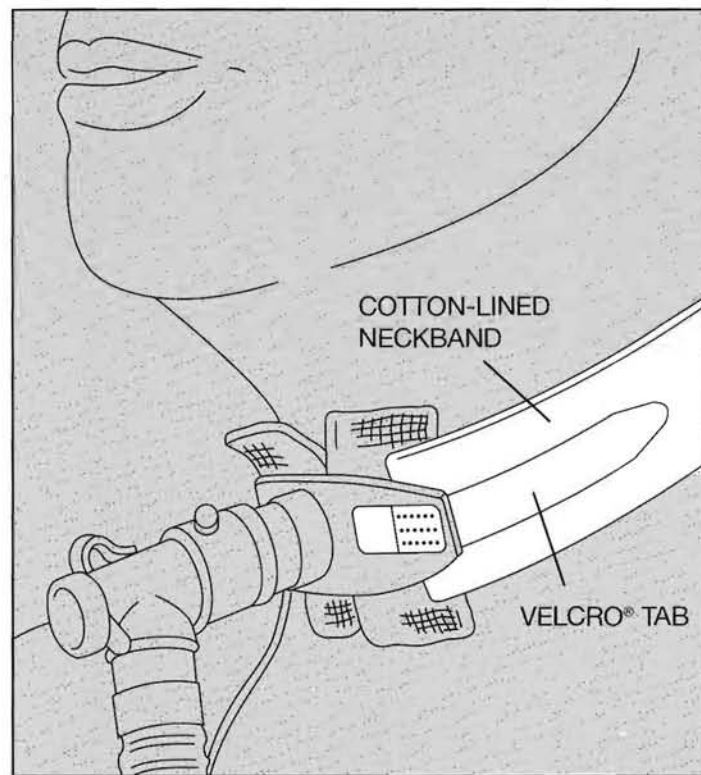

Ordering Information

DOUBLE COLLAR SET

PRODUCT NO. 43-1850

TRAKE-fit set includes two cotton-lined neckbands, (1-inch wide) with two wide (10mm) Velcro® hook tabs.

SINGLE COLLAR SET

PRODUCT NO. 43-1851

TRAKE-fit set includes one cotton-lined neckband, (1-inch wide) with two wide (10mm) Velcro® hook tabs.

NARROW TAB, DOUBLE COLLAR SET

PRODUCT NO. 43-1862

TRAKE-fit set includes two cotton-lined neckbands, (1-inch wide) with two narrow (6mm) Velcro® hook tabs.

PACKAGING

All products are packed non-sterile, in clear resealable plastic package. 25 per box, 100 per case.

*Velcro® is a Trademark of Velcro USA, Inc.

For maximum patient comfort

COMFIT™ Endotracheal Tube Holder

Significantly Decreases the Incidence of Self-Extubation*

- No tape on the ET tube; no tape on the patient
- May eliminate need for bite block
- Completely adjustable
- Wide cotton-lined neckband minimizes skin irritation, providing maximum patient comfort
- Minimal plastic loop around the ET tube allows access to the oral cavity
- Economical in two ways: low initial cost, no frequent changing
- Latex-free product
- Significantly decreases incidence of Self-Extubation*

Complies with ACLS Guidelines

ACLS Guidelines recommend the use of commercially available endotracheal tube holders as a substitute for tape.



The tapeless way to secure an endotracheal tube

Ordering Information

COMFIT

Product No. 42-1820

Plastic double securing loop, 1½" x 19" neckband.

Packaging

Non-sterile, 25 per box, 100 (4 boxes) per case.

* Reduction of Self-Extubation Using Comfit Tracheal Tube Holder
J. McAndrew RRT, Respiratory Care, Nov. 94, Vol. 39 No. 11.

ET Tube Fixation Made Simple

NEO-fit™ Neonatal Endotracheal Tube Grip

Provides quick, secure stabilization and allows for rapid release and readjustment of the ET tube

- Unique design eliminates the use of tape on the ET tube and minimizes the use of tape on the patient
- Bifurcated flex-pads provide secure facial contact
- Hook and loop strap allows rapid release and refit, making tube adjustment simple
- Fits 2.5mm to 4.0mm tubes, **not for use with 2.0mm or smaller tubes**



Ordering Information

NEO-fit

PRODUCT NO. 42-2540

Includes one plastic tube grip with adjustable hook and loop strap and adhesive foam foot pads, two foam "lollypop" securing tapes (optional use).



Packaging

Non-sterile, 20 per box, 80 per case (4 boxes).

Provides quick secure stabilization

Esophageal Balloon Catheter Set

For Intra-thoracic Pressure Monitoring

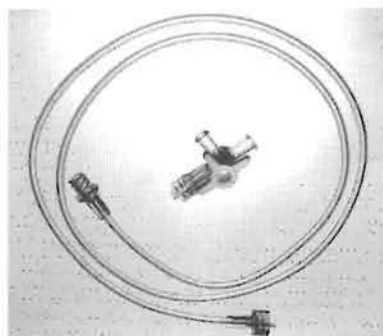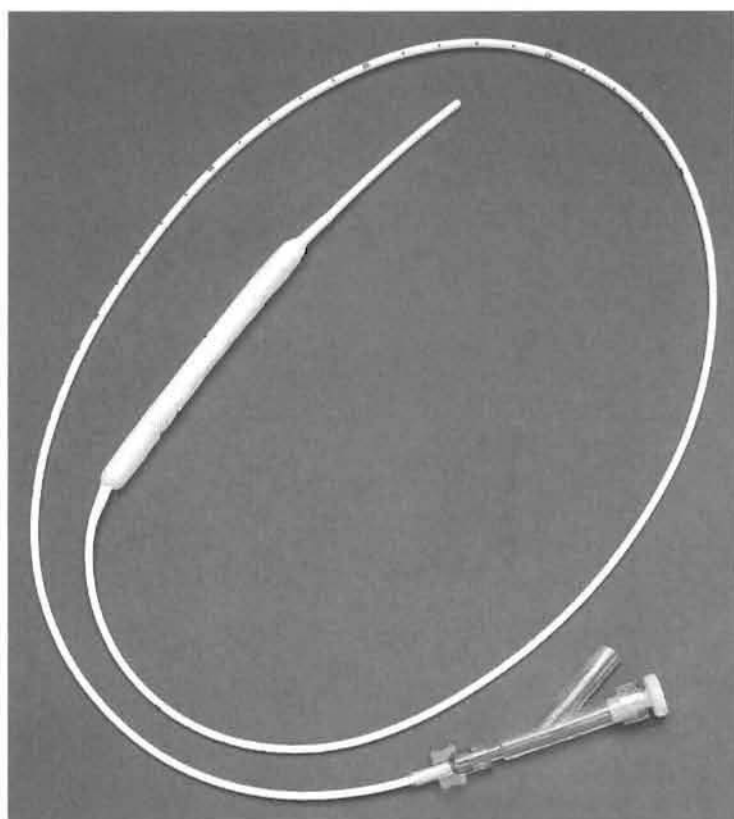

36" pressure monitoring extension and three-way stopcock

Ideal for monitoring ventilator and sleep apnea patients and calculating lung compliance.

- **Ease of insertion**

The catheter comes with a new easy-to-remove PTFE coated stylet to aid the insertion procedure

- **Comfort**

The 5 Fr. catheter is soft and pliable, offering maximum patient comfort over extended time periods

- **Ease of positioning**

The new proximal Y-connector allows for precise balloon positioning with the **stylet in place**

- **Depth markings**

The catheter is radio opaque with numerical depth markings to facilitate positioning

- **Ease of connection**

The catheter has a standard female luer connector

Ordering Information

Esophageal Balloon Catheter Set

Product No. 47-9005

Adult size. Includes graduated 5 Fr. catheter, 87cm. length, with soft latex-free balloon, 10cm length. PTFE coated stylet with rounded tip. 36" pressure monitoring extension, one three-way stopcock.

Packaging


Sterile, 5 sets per box, 20 sets per case.

INCA® Nasal Cannulae Set


The Clear Choice For Neonatal Nasal CPAP.

Five sizes for maximum stability and custom fit.

New NEO-prep Skin Barrier—Now Available

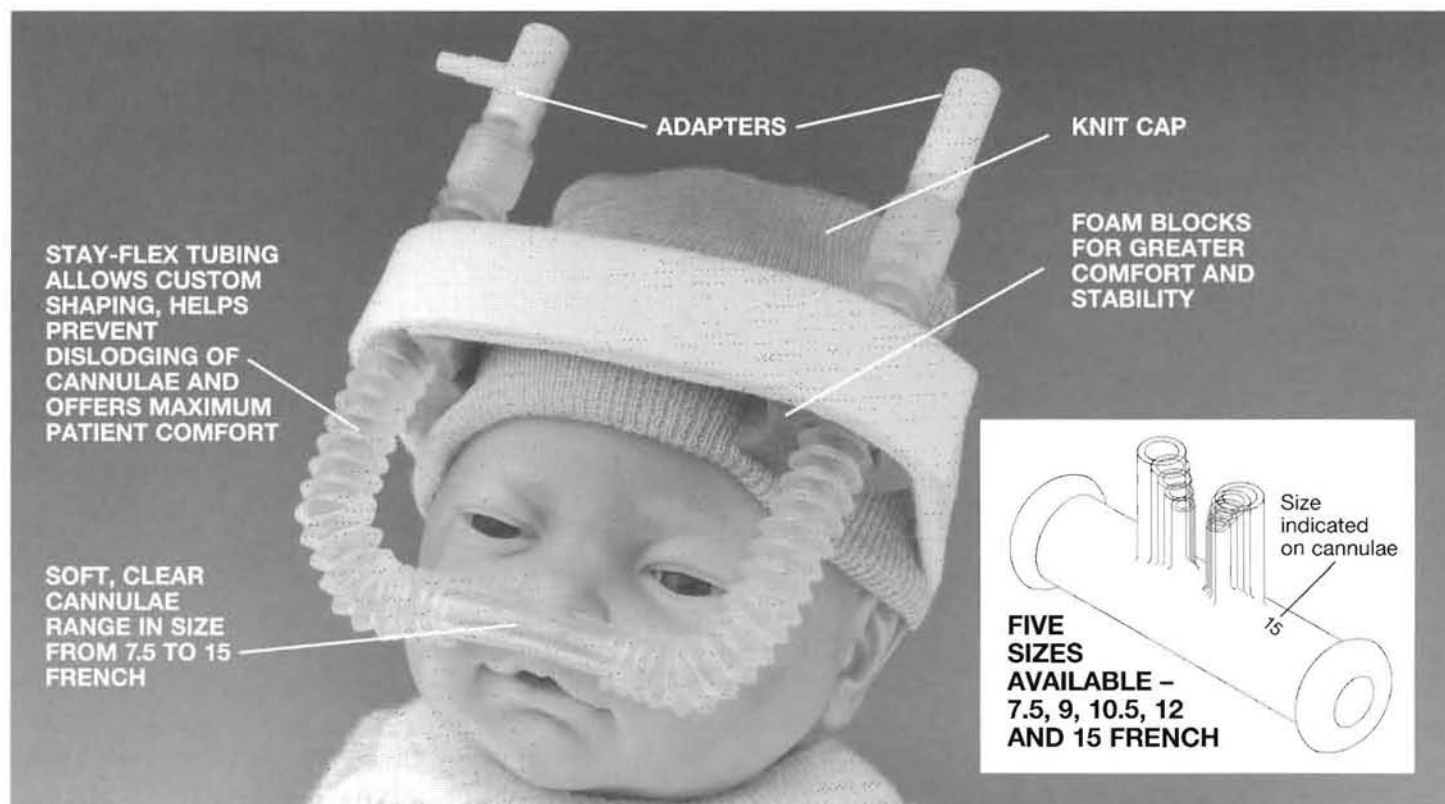

- Ultra-small 7.5 French size available
- Clear extra-soft cannulae
- Five close-increment sizes for more exact fit
- Stay-Flex tubing for custom shaping
- Foam block retainer for rapid placement and stability
- New NEO-prep minimizes septal breakdown ➤

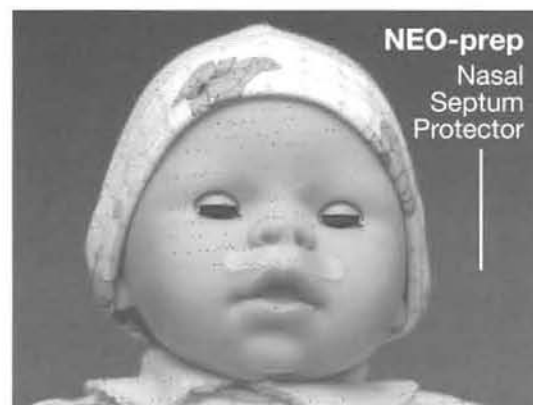

Ordering Information on Back Page

INCA® Nasal Cannulae Set Ordering Information

INCA® Complete Sets

PRODUCT NO. 44-2707

7.5 French nasal cannulae with connecting lengths of Stay-Flex tubing, knit cap with attached hook/loop band, 2 foam blocks, in-line pressure "T" adapter and 2 straight adapters. Complete set is packaged in easy to open tray. Nasal prongs packed sterile.

PRODUCT NO. 44-2709

Same as No. 44-2707 except with 9 Fr. cannulae.

PRODUCT NO. 44-2710


Same as No. 44-2707 except with 10.5 Fr. cannulae.

PRODUCT NO. 44-2712

Same as No. 44-2707 except with 12 Fr. cannulae.

PRODUCT NO. 44-2715

Same as No. 44-2707 except with 15 Fr. cannulae.

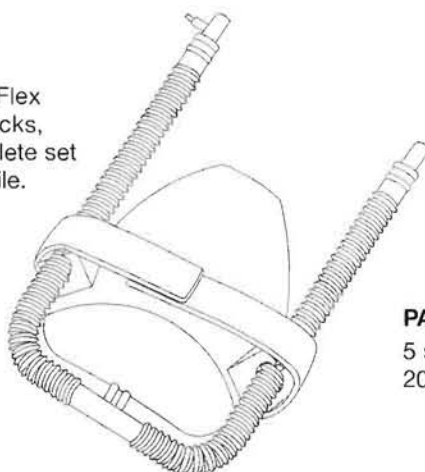

PACKAGING

5 sets per box,
20 (4 boxes) sets per case.

INCA® Replacement Sections

PRODUCT NO. 44-0707

7.5 French nasal cannulae with connecting lengths of Stay-Flex tubing, in-line pressure "T" adapter and 2 straight adapters. Set is packaged in easy to open tray. Nasal prongs packed sterile.

PRODUCT NO. 44-0709

Same as No. 44-0707 except with 9 Fr. cannulae.

PRODUCT NO. 44-0710


Same as No. 44-0707 except with 10.5 Fr. cannulae.

PRODUCT NO. 44-0712

Same as No. 44-0707 except with 12 Fr. cannulae.

PRODUCT NO. 44-0715

Same as No. 44-0707 except with 15 Fr. cannulae.

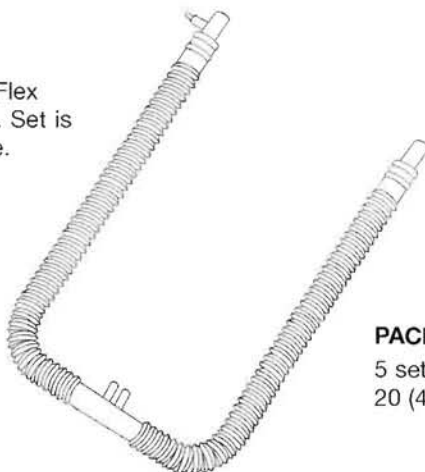

PACKAGING

5 sets per box,
20 (4 boxes) sets per case.

NEO-prep Skin Barrier

PRODUCT NO. 44-8000

Small pre-cut hypoallergenic skin barrier, non-sterile.

PRODUCT NO. 44-8001

Large pre-cut hypoallergenic skin barrier, non-sterile.

PACKAGING

8 per strip, 10 strips per box, 3 boxes (240) per case.



No. 44-8000

No. 44-8001

(Shown actual size)



CooperSurgical

95 Corporate Drive
Trumbull, CT 06611
Phone: 203-601-5200

International:
203-601-9818

*Keeping you at the forefront
of women's health care™*

IMA

Independent Medical Associates

11733 66th Street North
Suite 113
Largo, FL 33773

Phone: 888-548-4462
Fax: 888-548-1462
E-Mail: Info@I-MA.com