


H/S Catheter Set

For Hysterosonography and Hysterosalpingography

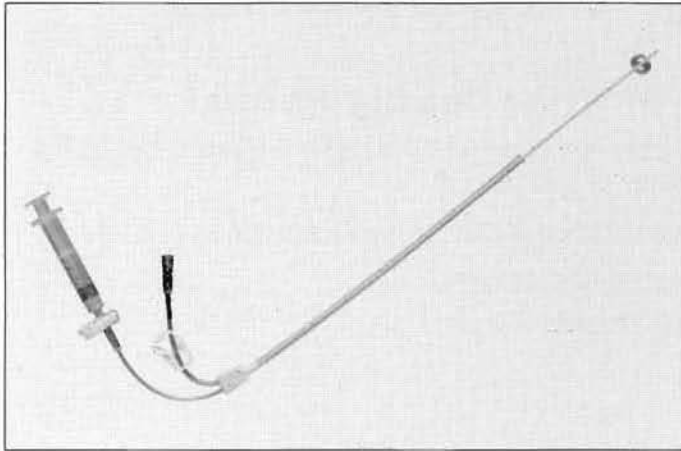

Features a unique placement sheath that eliminates the need for a tenaculum

PRODUCT NO. 61-5005

Flexible, radiopaque catheter, 5F, 30cm working length. Closed, rounded distal tip with side hole proximal to the tip. A non-latex balloon of 1.5cc capacity is mounted proximal to the side hole. The proximal end is bifurcated into one extension tube for balloon inflation and a second one for contrast or saline injection. Includes a 1.5cc controlled delivery syringe for proper balloon inflation. Intended for women who have not previously been pregnant.

PRODUCT NO. 61-5007

H/S Catheter Set, same as 61-5005 but 7F size with balloon capacity of 3cc. Includes a 3cc syringe. Intended for women who have previously been pregnant.

PACKAGING

Packed sterile in straight peel-pouch.
10 per box, 40 (4 boxes) per case.



H/S PROCEDURE TRAY

- Complete • Disposable • Sterile
- Cost Effective • Ideal for Office Use

PRODUCT NO. 61-5205

CONTAINS:
5 French H/S Catheter (61-5005) in sterile pouch
Fenestrated Drape • Disposable Speculum
Disposable Cervical Dilator • Two 3" x 3" Gauze Pads
Three Swab Sticks • One Wide Swab • Sterile Wrap
Inner Tray for Povidone - Iodine Solution
Povidone-Iodine Ointment • Lubricating Jelly Packet
Povidone-Iodine Solution (2 packets) • 36" Extension Tube
20cc Syringe with 18-Gauge Needle • Sanitary Napkin

PRODUCT NO. 61-5207

Same as 61-5205 except with 7F Catheter (61-5007).

PACKAGING

Packaged in sterile tray, 10 trays per case.

H/STYLET™


PRODUCT NO. 61-0020

Stainless steel stylet, 382mm long .020" diameter with luer lock hub. Used to provide additional catheter stiffness when indicated.

PACKAGING

Packed sterile in straight peel-pouch.
10 per box, 40 (4 boxes) per case.

Stainless Steel Instruments for Gynecologic Procedures

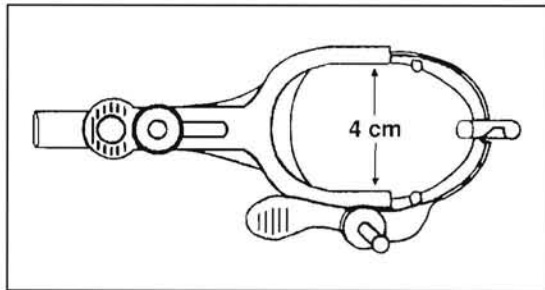


Vu-More Wide Rear-Opening Speculum

- Permits easy removal from the patient without disturbing the catheter or instrument
- The longer blade makes it suitable for obese patients

Product No. 63-902 Medium

Product No. 63-900 Large

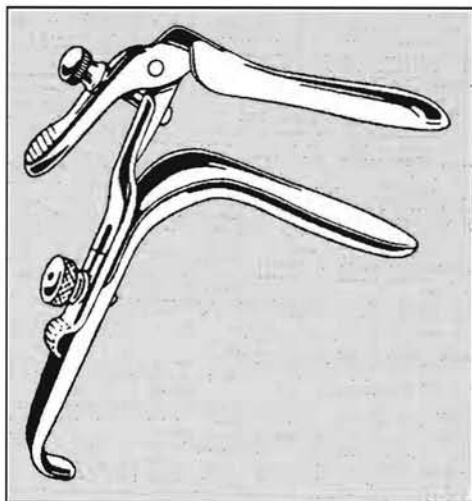

Weisman-Graves Open-Sided Speculum


- Speculum design permits easy removal without disturbing catheter or instruments
- Ideal time-saving device for Hysterosonography
- Easily adjustable to most anatomic needs
- Available with either left or right side opening

Product No. 391-667A Left side open – Medium

Product No. 391-668A Right side open – Medium

Product No. 391-667 Left side open – Large

Product No. 391-668 Right side open – Large

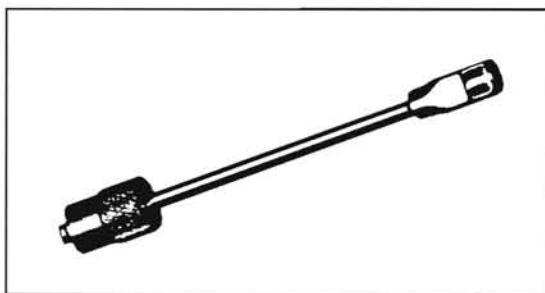

Standard Luer Needle Extension

- Allows use of less expensive standard size needles for paracervical blocks
- Malleable material can be bent to angle needed

Product No. 433-131A Reusable 3" extension

Product No. 433-132 Reusable 4" extension

Packaging

All of the products listed above are packaged individually in a clear plastic bag. Four per box. Non-sterile.

Distributed By:

IMA
Independent Medical Associates

Independent Medical Associates

11733 66th Street North
Suite 113
Largo, FL 33773

Phone: 888-548-4462
Fax: 888-548-1462
E-Mail: Info@I-MA.com