


Irrijet®

Provides fast, easy, forceful irrigation for removal of wound surface contaminants

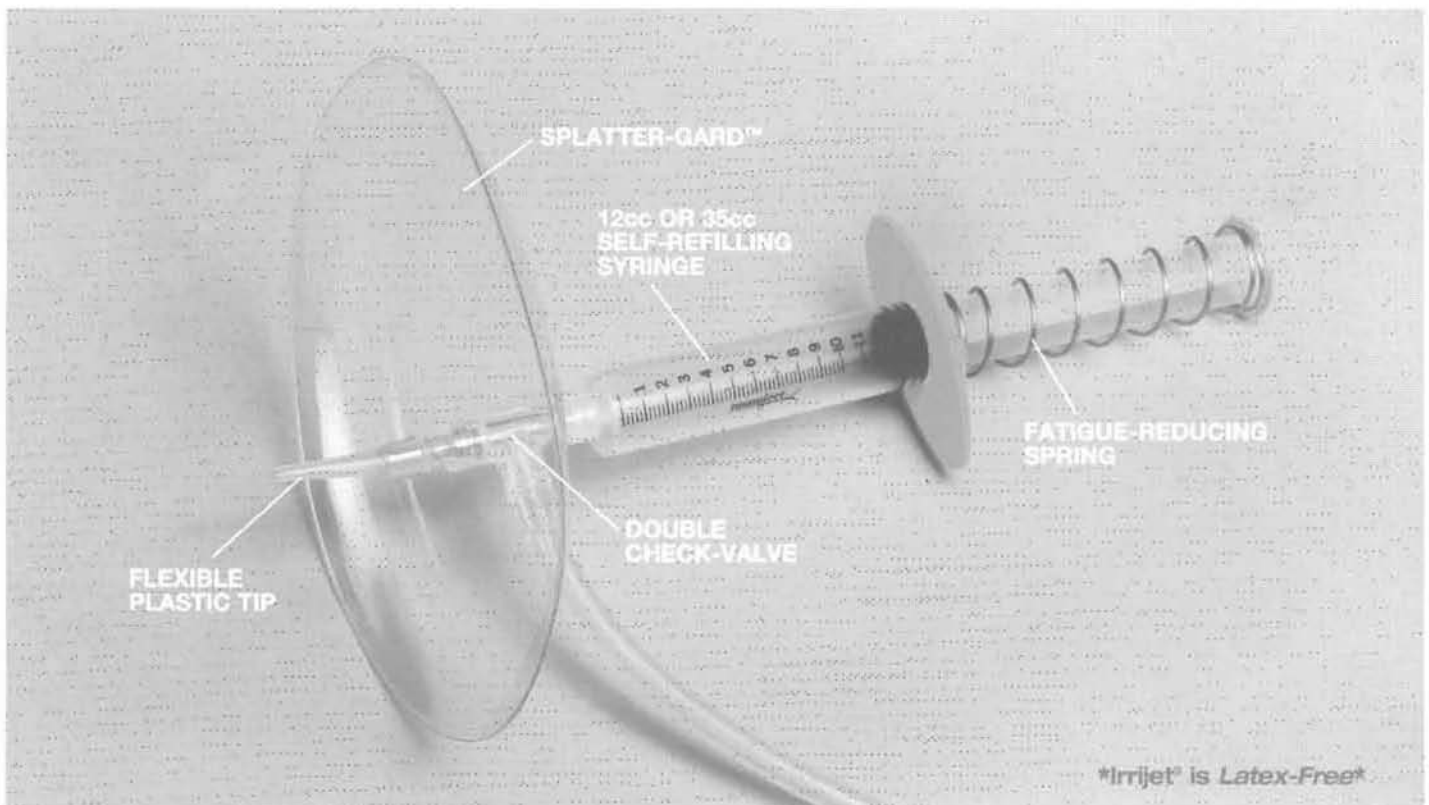

- **Splatter-Gard prevents contamination from body secretions—now with new concave design, provides maximum splash-back protection**
- **Soft, flexible plastic tip eliminates needle-stick for user, minimizes trauma to patient**
- **Irrigant is delivered with force equivalent to 8 psi**

The Irrijet is a unique, spring-loaded, self-refilling system that quickly and easily irrigates surface wounds.


Simple, one hand, thumb-controlled pressure provides an appropriately forceful stream (8 psi) directly to the desired area (or a gentler stream if needed). The soft flexible tip has an inside diameter equivalent to an 18-19 gauge needle.

The Irrijet is available in two sizes – 12cc or 35cc, and in single-use or multi-use models.

IRRIJET GP

PRODUCT NO. 55-0004

Spring-loaded syringe connected to a double check-valve, through a side arm and 6-ft. length of tubing to a spike for entry into a bag of irrigating solution. Double check-valve prevents retrograde flow. Alternate movement of the syringe plunger delivers irrigating solution as needed. Syringe capacity, 12 cc.

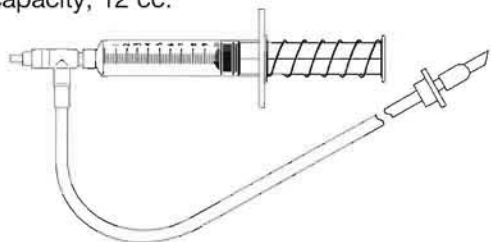


Packaging

Sterile, 5 per box, 20 per case (4 boxes).

Single-Use IRRIJET

PRODUCT NO. 55-2104

Same basic assembly as 55-0004 plus includes permanent flexible plastic jet tip and Splatter-Gard.

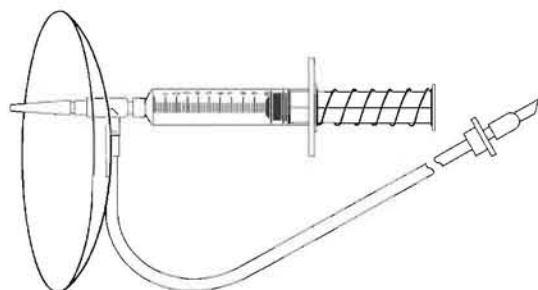

Packaging

Sterile, 5 per box, 20 per case (4 boxes).

Single-Use IRRIJET

PRODUCT NO. 55-2134

Same as 55-2104 but with 35cc capacity syringe.


Packaging

Sterile, 5 per box, 20 per case (4 boxes).

Single-Use IRRIJET Z

PRODUCT NO. 55-2004

Same basic assembly as 55-0004 plus one Zerowet Splashshield®.

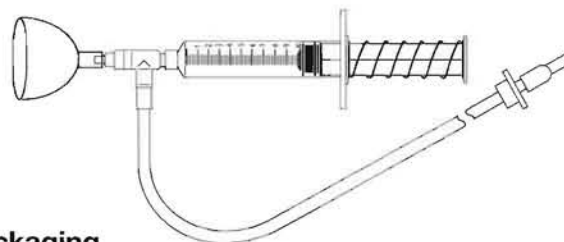


Packaging

Sterile, 5 per box, 20 per case (4 boxes).

Multi-Use IRRIJET

PRODUCT NO. 55-0504

Same basic assembly as 55-0004 plus one package of five individually sealed jet tips.

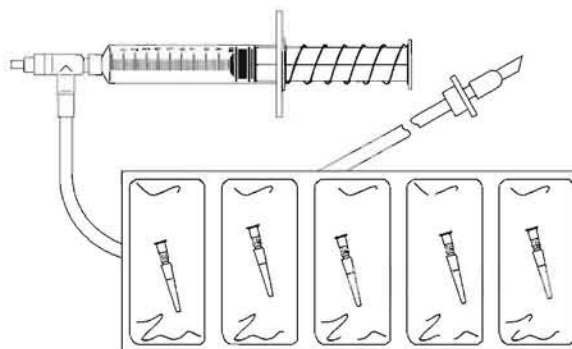

Packaging

Sterile, 5 per box, 20 per case (4 boxes).


Splatter Guard Tip and Cap Set

PRODUCT NO. 55-6100

Individually wrapped packages, each containing a Splatter-Gard, jet tip and cap.

Packaging

Sterile, 10 sets per box
4 boxes (40 sets) per case

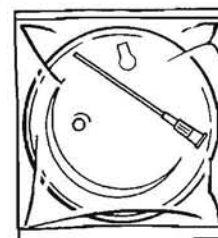

Distributed By:

IMA
Independent Medical Associates

Independent Medical Associates

11733 66th Street North
Suite 113
Largo, FL 33773

Phone: 888-548-4462
Fax: 888-548-1462
E-Mail: Info@I-MA.com