



INCA[®]

Infant Nasal CPAP Assembly

Your connection for Neonatal CPAP
and Bubble CPAP Systems

Features & Benefits

- Five anatomically proportional nasal cannulae sizes for a custom fit
- Clear, extra soft cannula helps maintain tissue integrity
- Stay-flex tubing provides custom shaping and maximum stability
- Foam block retainers provide secure placement and stability
- Chin strap accessory minimizes mouth leaks
- Available in both Complete and Replacement Sets
- Accessory Pack provides extra knit cap sizes, chin straps, hook fasteners and foam blocks


IMA
Independent Medical Associates



"Early bubble CPAP reduced delivery room intubations, (the number of) days on mechanical ventilation, postnatal steroid use and was associated with increased postnatal weight gain without any increased complications."

- Journal of Perinatology 2003 (23):195-199

Simple Bubble CPAP Setup and Components

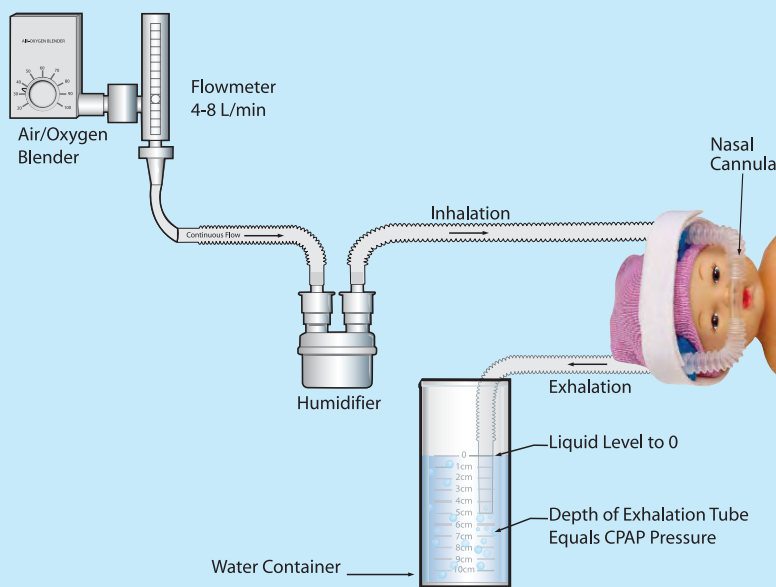

INCA[®] Infant Nasal CPAP Assembly, your connection for Neonatal CPAP (Continuous Positive Air Pressure) and Bubble CPAP Systems.



INCA Complete Set consists of a nasal cannula, two lengths of Stay-Flex tubing, two tubing adapters, one in-line "T" pressure adapter, a knit cap with hook/loop strap, two foam blocks and a hook/loop chin strap. One INCA Sizing Gauge is included in every box of five sets.

Complete Set P/N	Size
44-2707	7.5F
44-2709	9F
44-2710	10.5F
44-2712	12F
44-2715	15F



INCA Replacement Set includes the nasal cannula, two lengths of Stay-Flex tubing, two tubing adapters and one in-line pressure "T" adapter.

Replacement Set P/N	Size
44-0707	7.5F
44-0709	9F
44-0710	10.5F
44-0712	12F
44-0715	15F



INCA Accessory Pack includes five knit caps with hook/loop strap, five chin straps with ten hook fasteners (not attached) and ten foam blocks.


Accessory Pack P/N	Size
44-2001	X-Small
44-2002	Small
44-2003	Large
44-2004	X-Large



Nasal Sizing Gauge - For appropriate nasal prong fitting, use the INCA Sizing Gauge. Wipe gauge with sterile 70% isopropyl alcohol prior to each use. Hold the gauge up to the nasal area to determine appropriate cannula size.

Nasal Sizing Gauge P/N
025-001

**NOTE: Do not place sizing gauge in neonate's nares.*

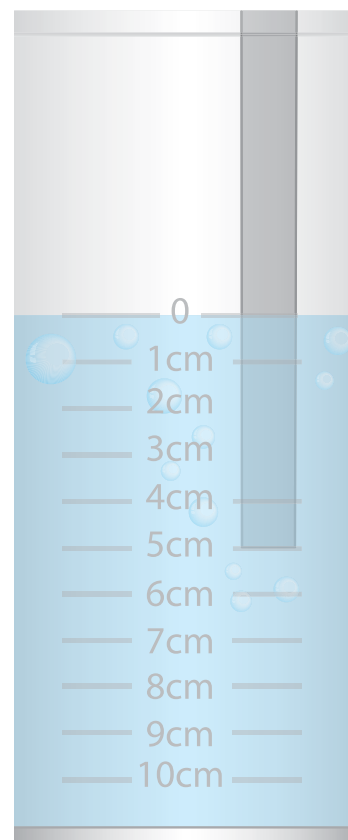

IMA
Independent Medical Associates

11733 66th St. N #113 Largo, FL 33771 • 727.548.4462 • 888.548.4462 • Fax 888.548.1462
www.I-MA.com • ©2007 CooperSurgical, Inc.

Form # 81256 06/07

INCA[®] is a registered trademark of CSI.