

# The Most Advanced Method to Diagnose and Treat Patients with Chronic Low Back Pain

**VUECATH<sup>®</sup>**  
*System*

Spinal Endoscopic System

Accurate Diagnosis,  
Effective Treatment





Dual Lumen Catheter  
now with HALOSCOPE<sup>®</sup>  
fiberoptics



# VueCath® Spinal Endoscopic System

## Vivid Image from a Reusable Fiberscope

The VueCath Spinal Endoscopic System makes it possible for the physician to directly visualize the lumbar and sacral spine anatomy under local anesthesia. Used diagnostically, the System allows the physician to view the pathology and make an accurate diagnosis; used by the pain management physician, it allows for effective treatment.

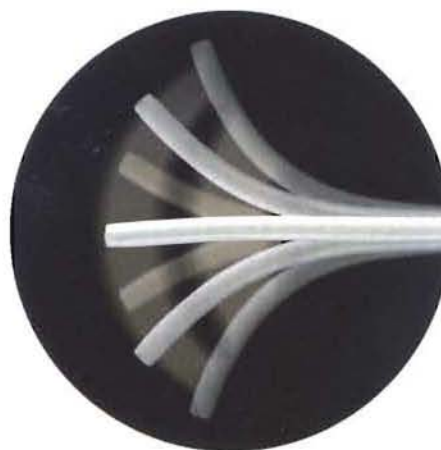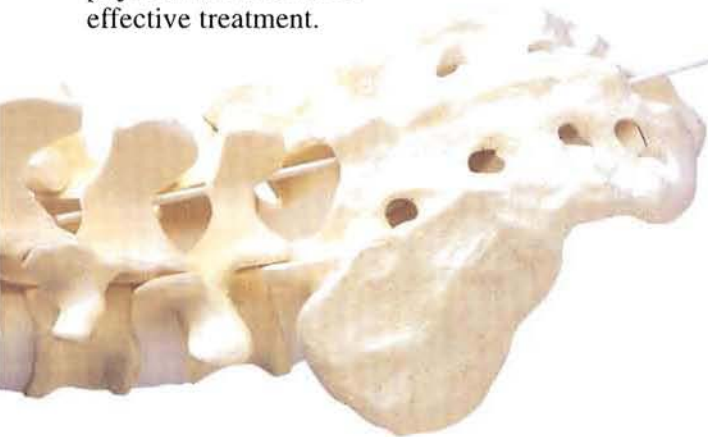

## Accurate Diagnosis, Effective Treatment

- Small, flexible, reusable 1.0mm fiberscope has 10,000 pixels, producing a **vivid picture**;
- Single-use catheter **eliminates cross contamination** concerns, and eliminates need to clean catheter;
- Ergonomically designed **steerable catheter** allows precise tip placement;
- Small diameter of catheter allows for **easy entry into epidural space**;
- Standard camera coupler **enables scope to be adapted to most existing video systems**;
- Round, **atraumatic**, distal tip allows visualization while the scope is withdrawn into catheter;
- Fluid port provides secure fitting **for approved drug delivery** to the site of the pain, preventing leakage;
- Auxiliary port **offers an additional lumen for administration of fluids**.
- A highly effective tool for **mechanical blunt dissection** in the epidural space

## Ordering Information

Cat. No.	Description
S-004	VueCath HALOSCOPE® Fiberscope - 1.0mm OD
S-107	VueCath Fiberscope Coupler - 7X with C-mount.
S-111	VueCath Eye Piece Adapter with C - mount, to be used with cat # S-107
S-113	VueCath V- mount Adapter, to be used with cat # S-107
C-005	VueCath Steerable Catheter, Sterile, Contains 5 Catheters
C-010	VueCath Steerable Catheter, without Detents, Sterile, Contains 5 Catheters
Y-001	VueCath Endoscopic Disposable Accessory Kit, Sterile, Contains 5 kits
Y-005 CE	VueCath Endoscopic Disposable Accessory Kit (International CE mark use) Sterile, Contains 5 kits

For more information, contact EBI today or call **1-800-526-2579**.

For a current list of VueCath System products, contact your representative.



11733 66th Street North  
Suite 113  
Largo, FL 33773

Distributed by:  
Independent Medical Associates

Phone: 888-548-4462  
Fax: 888-548-1462  
E-Mail: Info@I-MA.com

This product and its methods for use are covered under one or more of the following U.S. Patents: Des. 343,678; Des. 349,340; 5,342,299; 5,399,164; 5,423,311; 5,437,636; 5,531,687; 5,624,397; 6,010,493; 6,464,682; 6,470,209. Other U.S. and foreign patent rights pending and/or issued.

### CAUTION:

Federal law (USA) restricts this device to sale by or on the order of a physician. EBI is a registered trademark of Electro-Biology, Inc. VueCath and HALOSCOPE are registered trademarks of, and are distributed by, EBI, L.P. For full prescribing information, contact EBI, L.P., a subsidiary of Biomet.

3013 5M ©2003 EBI, L.P. All rights reserved.



EBI, L.P.  
100 Interpace Parkway  
Parsippany, NJ 07054  
www.ebimedical.com  
1-800-526-2579

A BIOMET COMPANY