

Soft  
&  
Strong

# Feth-R-Kath™



**Stingray™ Connector**  
NEW Stingray™ Connector is packaged with all Epimed Catheters.

Made in the USA by Epimed International

## Benefits:

- Kink-Resistance
- Soft, Atraumatic, & Flexible Tip
- Welded Safety Ribbon
- Radio-Opacity
- Natural Lubricity of Nylon
- Packaged with Stingray Connector™


## Suggested Introducers:

- Tuohy      ▶ RX™
- Crawford    ▶ RK™

19g x 24 inches Re-order#: 157-1924

## Specifications:

The 19g Feth-R-Kath™ comes in a convenient 24" (61cm) length for ease of handling and enhanced direct-ability. This catheter features a proprietary nylon blend coating for ideal insertion and removal characteristics as well as shear resistance. Feth-R-Kath™ can be manipulated with confidence into the epidural space with minimal chance of catheter problems.

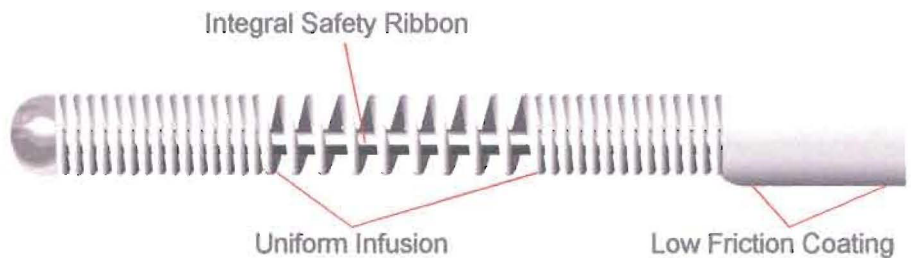

SPIROL®

Soft & Strong

# Spirol® Catheter

Spring Reinforced Epidural Catheter



### Stingray Connector™

NEW Stingray Connector™ is packaged with all Epimed Catheters.

Made in the U.S.A. by Epimed International

### Features:

The new SPIROL® catheter provides a confident yet gentle means of administering medication. Increased softness and flexibility ensures a more comfortable experience for patients. The SPIROL® is the perfect balance of strength and flexibility in delivering medication.

### Benefits:

- ▶ Kink-Resistance
- ▶ Natural Lubricity of Nylon
- ▶ Soft, Atraumatic, & Flexible Distal Tip
- ▶ Columnar Strength
- ▶ Superior Flow Rates
- ▶ Radio-Opacity
- ▶ Flashback Windows
- ▶ Thermal Formed Tip


### Specifications:

Traditional epidural catheters offer compromises - strength means stiffness yet softness risks kinking and breakage. The SPIROL® catheter offers the ideal balance of strength and softness. The spring coil is encapsulated in a proprietary nylon blend coating. The natural lubricity of nylon slides through tissue as opposed to urethane coated spring catheters which tend to grab during insertion or removal. From the soft flexible tip, to the firm kink resistant body, the SPIROL® is designed to deliver a new measure of safety and fault-free performance.

*For use with a 17g Epidural Needle.*


*20g version and various kit configurations available soon.*

### Catheter Tips

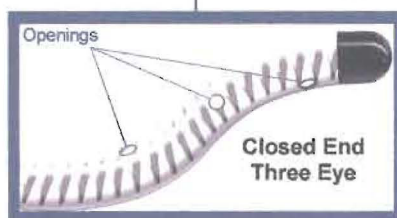

19g x 36 inch Closed End 3 Eye

Re-order#: 166-1936

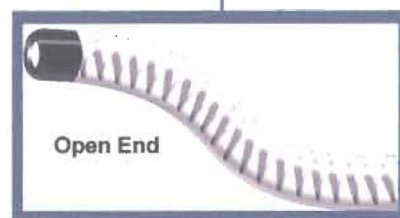

19g x 36 inch Open End

Re-order#: 165-1936