

# PRECISA M

INNOVATIVE PRODUCTS FOR BREAST BIOPSY

Thanks to the needle removable outer cannula, it is possible to perform either multiple histological biopsies or a double micro-histological and cytological biopsy.

The new 13G size is added to the traditional gauges.

#### VERSATILITY:

It is designed to perfectly fit in with the most recent imaging diagnostic instruments (Stereotaxis - prone biopsy table - Echographers) through particular brackets (on request).

Both colouring liquids use or a proper metal-guidewire positioning can be performed with this needle during the histological biopsy.

#### PRECISION and SAFETY:

1cm and 2cm pre-set penetrations for specific histological biopsies.

Division of different biopsy stages and control of advancement are performed safely even at different speed thanks to the "SAFETY BLOCK".

- Echogenic marker
- Depth marks
- Extra stylet with pyramidal tip
- A two rings grip
- Stereotaxis spacer


CODE	Diameter	Length	Box
PRE1310	13G	100mm	5
PRE1410	14G	100mm	5
PRE1415	14G	150mm	5
PRE1420	14G	200mm	5
PRE1610	16G	100mm	5
PRE1810	18G	100mm	5

#### BRACKETS

CODE	Description	Box
ADFPRE	Bracket for Fischer table	1
ADPRE	Bracket for Larad table	1
ADGPRE	Bracket for Giotto* table	1

\*Giotto is an I.M.S. srl's trade-mark

Patent Pending

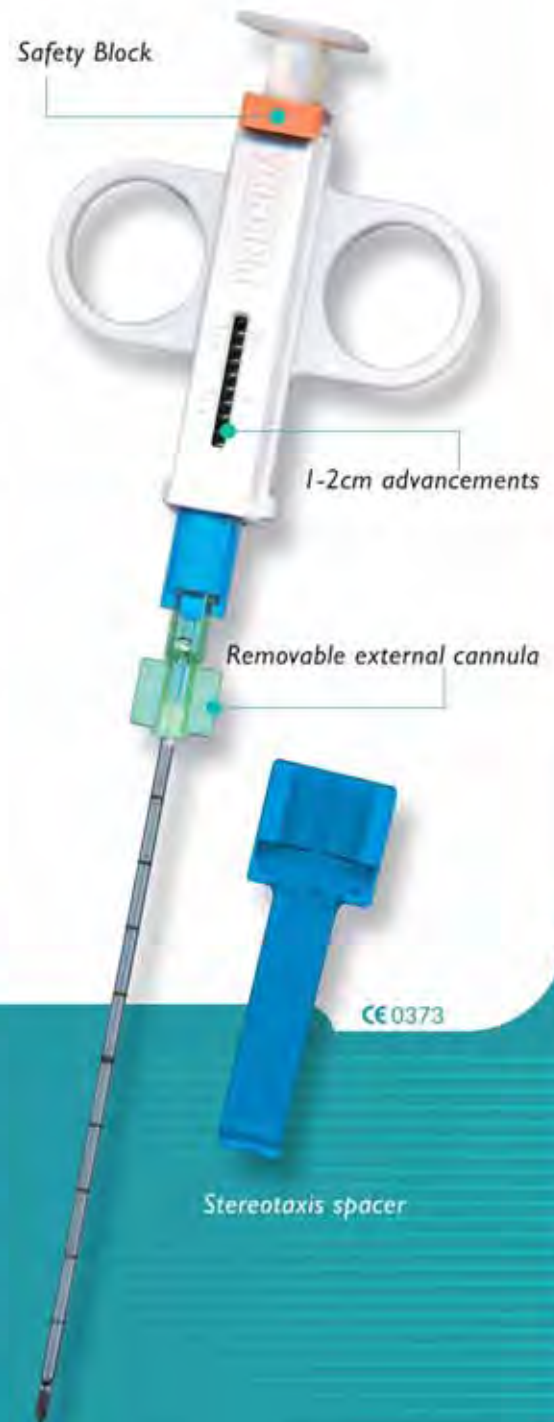

# IMA

Independent Medical Associates

www.I-MA.com  
info@I-MA.com  
p: 888.548.4462  
f: 888-548-1462