

PROPPER

The difference

INTEGRATED

you can

FIBER OPTIC

see...with.

LARYNGOSCOPES

propper



Brilliant.

IN MORE WAYS THAN ONE.

Propper Integrated Fiber Optic Laryngoscopes provide you with one brilliant advantage after another. Brighter light source, patient comfort, sterilization, flexibility and ease of cleaning.

With over 60 years of experience and service to the medical profession, Propper has a pioneering history of providing the most innovative sterilization products on the market, backed by our commitment to quality, reliability and service.

Take a good look at Propper Integrated Fiber Optic Laryngoscopes, and see for yourself why more and more professionals are choosing Propper.

*Patented
Switch
Mechanism
in Handle.*



*High Intensity
Reflector
Bulb
in Handle.*

BRILLIANT LIGHT:

High-intensity reflector bulbs provide you with a whiter, brighter light. A yellow tinge to the light alerts you to the need to insert fresh batteries.

BRILLIANT DESIGN:

We put the bulb in the handle, not the blade! And that helps with patient comfort because the blade does not heat up like other units where the bulb is in the blade. What's more, there are no wires to repair because there is no electrical conductivity between the blade and handle. Blades simply snap into place activating Propper's patented switching mechanism* - a mechanism which insures reliable long term performance due to its positive contact design.


Propper Laryngoscopes are easy to clean, because the fiber optic bundle is a permanent part of the stainless steel blade. There is no need to take the blade apart for cleaning like many other blades. **And, as a leader in the sterilization field, both the blade and outer compartment of the Propper handle are manufactured to be sterilized with either STERRAD® or STERIS® sterilization systems**.** All Propper Laryngoscopes can be sterilized with any approved EO gas sterilization cycle as well.

Propper Integrated Fiber Optic Laryngoscopes come in three handle size variations with numerous blade types and sizes in each set.

When you choose a Laryngoscope, make sure you look for the differences you can see, and also the ones you can't. Like Propper engineering excellence, reliability and quality, and a sales, service and distribution network that is dedicated to meeting all of your needs.

* U.S. Patent number 5,542,905

** STERRAD is a registered trademark in the United States of Johnson & Johnson. STERIS is a registered trademark in the United States of Steris Corporation.

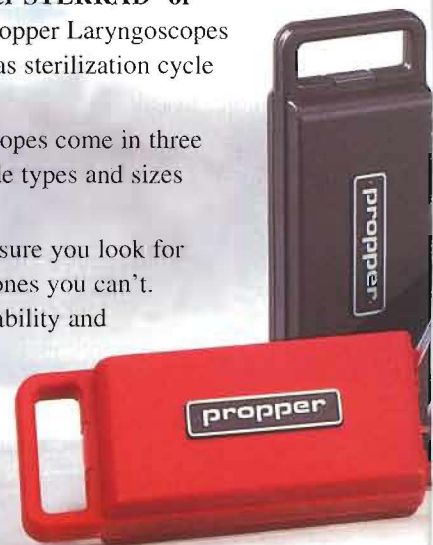

PROPPER OFFERS A COMPLETE LINE OF FIBER OPTIC LARYNGOSCOPES

COMPLETE SET 199177

Contains 5 Macintosh blades, 5 Miller blades, 3 Wisconsin blades, 2 Propper blades, and 3 handles in rigid black attache case.

LARGE SET 199176





Contains 5 Macintosh blades, 5 Miller blades, 2 Propper blades and 2 handles in rigid black attache case.

SMALL SET 199175

Contains 1 each of Propper #1, Macintosh #3, and Miller #3 blades, plus a large handle and a small handle in a compact soft case.

**EMERGENCY SETS
19917502/19917503**

Contains 1 each of Propper #1, Macintosh #3, and Miller #3 blades, plus a large handle and a small handle in a rigid red or black compact case.

Macintosh Blades		
	Size	Reorder No.
	1	199161
	2	199144
	3	199145
	4	199146
	5	199162
Miller Blades		
	Size	Reorder No.
	0	199156
	1	199157
	2	199158
	3	199159
	4	199160
Wisconsin Blades		
	Size	Reorder No.
	1	199163
	2	199164
	3	199165
Propper Blades		
	Size	Reorder No.
	0	199153
	1	199147



One

BULB IS ALL YOU NEED.

Because the bulb is in the handle, only one size bulb is needed to light any Propper fiber optic laryngoscope blade.

One bulb (reorder No. 199407) fits all three 2.5V handles. Unscrews for easy replacement.



Propper Fiber Optic Laryngoscope handles for fiber optic blades.
 No. 199150 uses 2 "C" batteries
 No. 199152 uses 2 "AA" batteries
 No. 199578 uses 2 "AA" batteries

IMA

Independent Medical Associates

Distributed by:

Independent Medical Associates

11733 66th Street North
 Suite 113
 Largo, FL 33773

Phone: 888-548-4462
 Fax: 888-548-1462
 E-Mail: Info@I-MA.com

ORDERING INFORMATION

Available through local Anesthesia dealers

BLADES	DIMENSIONS		REORDER NUMBER
	Overall Length (mm)	Vertical Opening (mm)	
Macintosh			
1	93.0	15.0	199161
2	115.0	16.0	199144
3	135.0	16.5	199145
4	160.0	17.0	199146
5	176.0	18.0	199162
Miller			
0	78.0	9.4	199156
1	102.0	10.4	199157
2	152.0	11.4	199158
3	192.0	12.8	199159
4	207.0	15.2	199160
Wisconsin			
1	110.0	14.6	199163
2	130.0	16.8	199164
3	160.0	17.5	199165
Propper (straight blades for children with large tongues)			
0	83.0	9.3	199153
1	95.0	9.3	199147
Sets			
Small Emergency Set (Hard Case-Red)			19917502
Small Emergency Set (Hard Case-Black)			19917503
Small Set (Soft Zippered Case)			199175
Large Set			199176
Complete Set			199177
Handles (batteries included)			
2.5V — Uses 2 "C" batteries			199150
2.5V — Uses 2 "AA" batteries			199152
2.5V — Uses 2 "AA" batteries — Stubby			199578
Bulb			
2.5V — Reflector, high-intensity			199407
Cases			
Small Emergency Set (Hard Case-Red)			19917902
Small Emergency Set (Hard Case-Black)			19917903
Small Set (Soft Zippered Case)			199179
Large & Complete Sets Attache Case			199180

In addition to receiving help from your Propper field representative, you can also use the Sterilization Hotline (800) 221-2280. Propper scientists are available to answer your questions from 9 a.m. to 5 p.m., EST.

A Schuman Company

propper

propper manufacturing company inc.

36-04 Skillman Avenue
 Long Island City, NY 11101
 Customer Service: (800) 832-4300
 www.proppermfg.com

©1997 Propper Printed in USA 01/00